

TS-150-S Series Active Vibration Isolation



TS 150-S

Actively damped frequency range: **0.7 – 1000 Hz**

Load capacity: **0 - 150 kg**

Dimensions: **400 x 450 x 75 mm**

This dynamic vibration isolation system achieves, in a very small volume, better isolation than is possible with most any isolation system. This new technology provides isolation from building vibrations, and also from vibration sources placed on the tabletop. A precise automatic level adjustment mechanism allows convenient handling of different loads without any user intervention.



The TS-150 S offers ultimate performance for a wide range wide range of high resolution instruments. e.g. *SPM, interferometer, microscopy*

- Isolation begins at 0.7 Hz, increasing rapidly to 40 dB at 10Hz
- No low frequency resonance
- 500 times stiffer than air tables.
- Isolation of all six possible translational and rotational modes of vibration.



HERZ CO.,LTD.

18/F., YCS Bldg., 5-1, Sakae-Cho, Kanagawa-Ku,
Yokohama, Kanagawa 221-0052 Japan
Phone : +81 (0 in JP) 45-450-2211 Fax : +81 (0 in JP) 45-450-2221
e-mail : sales@herz-f.co.jp <http://www.herz-f.co.jp>




HERZAN L.L.C.

23151 Alcalde Dr. Unit B-3 Laguna Hills CA 92653 U.S.A
Phone : 949-363-2905 Fax : 949-363-2906
e-mail : sales@herzan.com <http://www.henzan.com>

Transmissibility

TS-150-S

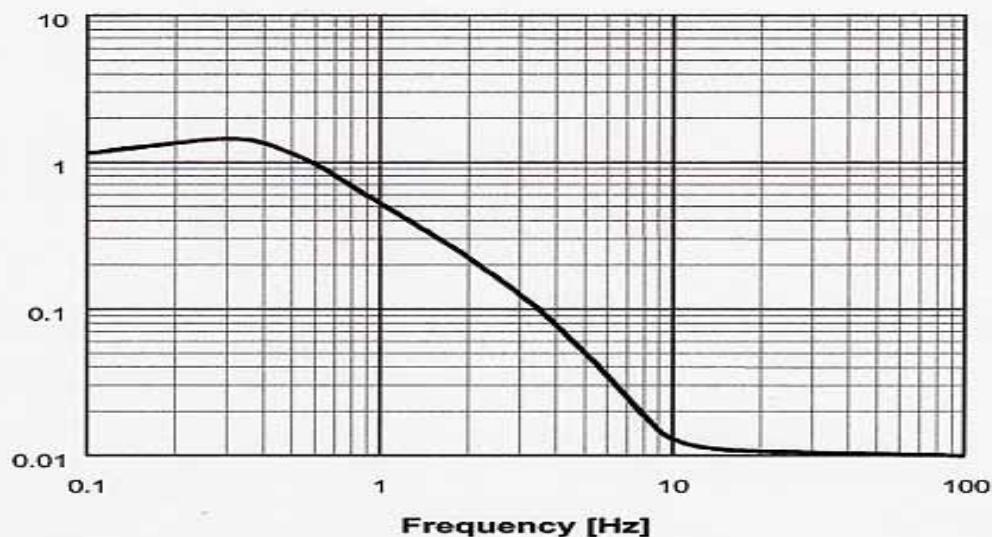

Isolation: dynamic 0.7 Hz to 1000 Hz, passive beyond 1000 Hz.

Transmissibility: Above 10 Hz transmissibility < 0.01 (-40dB)

System Noise: less than 50nG/ $\sqrt{\text{Hz}}$ from 0.1-200 Hz in any direction

Maximum Load: 150 Kg

Size: 400 x 450 x 75 mm

Weight: 18 kg

ELECTRICAL

Safety class: 1

Power Consumption: typically 10 W, maximum 50 W (70 VA)

Input Voltage: 115-230V AC $\pm 10\%$, 50-60Hz

Monitor Signal: A multiplexed signal for display on oscilloscope shows vibration levels with and without isolation. A signal of 1mV corresponds approximately to a velocity of 0.5 μ m/sec at the table top.

TABLE TOP

Size: 400 x 450 mm

Material: Aluminium

Also available as standard ferromagnetic plate

TS-150 : Examples of application

